



FavorPrep™ Plasmid Extraction Maxi Plus Kit

For Research Use Only

Kit Contents:

	FAPMX 002 (2 preps)	FAPMX 010 (10 preps)
PEQ Buffer	30 ml	135 ml
PM1 Buffer	42 ml	215 ml
PM2 Buffer	42 ml	215 ml
PM3 Buffer	42 ml	215 ml
PW Buffer	65 ml	270 ml + 60 ml
PEL Buffer	32 ml	215 ml
RNase A Solution	100 µl	480 µl
PMX Column	2 pcs	10 pcs

* Preparation of PM1 Buffer for the first use:

Cat. No:	FAPMX 002 (2 preps)	FAPMX 010 (10 preps)
Volume of RNase A Solution for PM1 Buffer	84 µl	430 µl

Specification:

Technology: Anion-exchange chromatography (gravity-flow column)

Lysate clarification: Centrifugation

Sample Size: 120~240 ml of bacteria for high-copy number or low-copy number plasmid

Plasmid or constructs range: 3 kbp~150 kbp

Binding Capacity: 1.5 mg/Maxi Column

Important Notes:

1. Store RNase A Solution at -20°C upon receipt of kit.
2. Add indicated volume of RNase A Solution into PM1 buffer, mix well and store the PM1 buffer at 4°C.
3. If precipitates have formed in PM2 Buffer, warm the buffer in 37°C water bath to dissolve precipitates.
4. Pre-chill PM3 Buffer at 4°C before starting.

Additional Requirements:

1. 50 ml tubes
2. Refrigerated centrifuge capable of $\geq 5,000$ xg and the centrifuge tube suitable for the centrifuge rotor
3. Isopropanol
4. 70% ethanol
5. TE buffer or ddH₂O

General Protocol:

Please Read Important Notes Before Starting Following Steps.

Harvest bacterial cells


1. Harvest the cells by centrifugation at 4,500~6,000 xg at 4°C for 10 mins and discard the supernatant.

Equilibrate PMX Column

2. Hold a PMX Column above a waste tank.
3. Equilibrate the PMX Column by applying 10 ml of PEQ Buffer. Allow the column to empty by gravity flow and discard the filtrate.

Cell lysis and lysate neutralization

4. Add 16 ml of PM1 Buffer (RNase A added) to resuspend the cell pellet by vortexing or pipetting.
5. Add 16 ml of PM2 Buffer and mix gently by inverting the tube 5 times.
-Do not vortex to avoid shearing genomic DNA.
6. Incubate the sample mixture for 5 mins at room temperature until lysate clears.

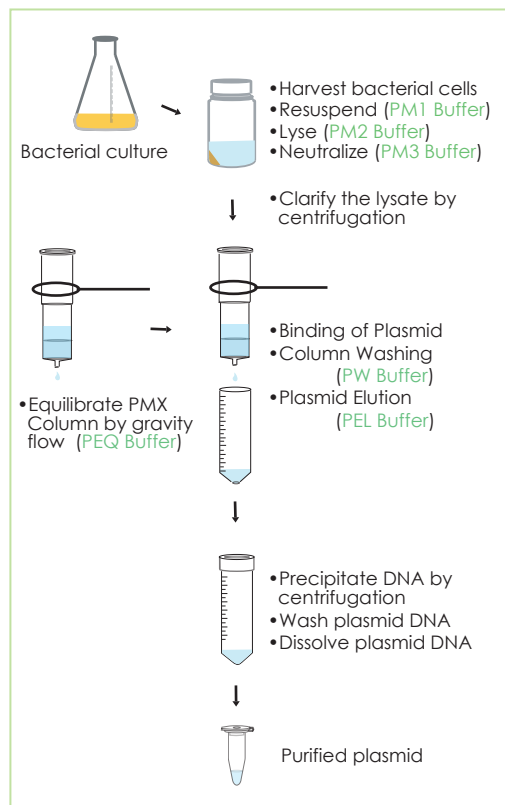

7. Add 16 ml of chilled PM3 Buffer and mix immediately by inverting the tube 10~15 times to neutralize the lysate. (Do not vortex!)

- Note:
- Make sure the density of cultured cell is optimal, the buffers volume (PM1, PM2, PM3) should be increased proportionally to the culture volume.
(ex. culture volume, 120~240 ml: PM1, 16 ml; PM2, 16 ml; PM3, 16 ml
culture volume, 240~480 ml: PM1, 32 ml; PM2, 32 ml; PM3, 32 ml)
 - Make sure cell pellet be suspended completely within Buffer PM1.
 - Mix the sample mixture completely after adding Buffer PM2 and Buffer PM3.

Clarify lysate by centrifugation

8. Centrifuge the tube at $\geq 5,000$ xg at 4°C for 20 mins. (Preferably centrifuge the tube at 15,000~20,000 xg at 4°C for 15 mins.)

- If the supernatant still contains suspended matter, transfer the supernatant to a clean centrifuge tube and repeat this centrifugation step.

Binding of plasmid

9. Transfer the supernatant from step 8 to the equilibrated PMX Column. Allow it to flow through the PMX Column by gravity flow and discard the filtrate.

Wash PMX Column

10. Wash the PMX Column by applying 30 ml of PW Buffer. Allow PW Buffer to flow through the PMX Column by gravity flow and discard the filtrate.

Elution

11. Hold the PMX Column above a clean 50 ml centrifuge tube (not provided). Add 15 ml of PEL Buffer to the PMX Column to elute the plasmid by gravity flow.

Precipitate plasmid DNA

12. Add 0.75X volume of room temperature isopropanol to the eluate and mix well by inverting the tube 10 times. (ex: add 11.25 ml isopropanol to 15 ml eluate)

- Note! Make sure that isopropanol be mixed thoroughly with eluate before centrifugation.

13. Centrifuge the tube at $\geq 5,000$ xg at 4°C for 30 mins. (Preferably centrifuge the tube at 15,000~20,000 xg at 4°C for 20 mins.)

Wash and dissolve plasmid DNA

14. Carefully remove the supernatant and wash the plasmid pellet with 5 ml of room temperature 70% ethanol.

15. Centrifuge the tube at $\geq 5,000$ xg at 4°C for 10 mins.

16. Carefully remove the supernatant and invert the tube on paper towel for 3 mins to remove residual ethanol. Air-dry the plasmid pellet until the tube is completely dry. (Or incubate the plasmid pellet at 70°C for 10 mins.)

17. Dissolve the plasmid pellet in a suitable volume (≥ 300 μ l) of TE or ddH₂O.

- Note!
- Do not lose the DNA pellet when discard the supernatant.
 - Make sure the DNA pellet adhesive lightly on the centrifuge tube.
 - If the DNA pellet loose from tube, repeat the precipitation step again.
 - Make sure the DNA is dissolved completely before measure the concentration.

Troubleshooting

Low yield

Bacterial cells were not lysed completely

- Too many bacterial cells were used.
- After PM3 Buffer addition, break up the precipitate by inverting.
- DNA failed to precipitate or DNA pellet was lost after precipitation.
- DNA pellet was insufficiently redissolved.

Purified DNA dose not perform well in downstream applications

RNA contamination

- Make sure that RNase A had been added in PM1 Buffer at the first use.
- RNase A is not properly preserved.
- Too many bacterial cells were used, reduce the sample volume.

Genomic DNA contamination

- Do not use overgrown bacterial culture.
- During PM2 and PM3 Buffer addition, mix gently to prevent genomic DNA shearing.
- Lysis time was too long (over 5 mins).

Too much salt residual in DNA pellet

- Wash the DNA pellet twice with 70% ethanol.