

## Kit Contents:

Cat. No. (preps)	FABGK000-Maxi (2 preps)	FABGK003 (10 preps)	FABGK003-1 (24 preps)
Proteinase K (Liquid)	1050 µl × 2	10.5 ml	13 ml × 2
FABG Buffer	22 ml	110 ml	265 ml
W1 Buffer* (Concentrate)	6.5 ml	33 ml	88 ml
Wash Buffer** (Concentrate)	3 ml	20 ml	40 ml
Elution Buffer	6 ml	30 ml	60 ml
FABG Maxi Column	2 pcs	10 pcs	24 pcs
Elution Tube (50 ml tube)	2 pcs	10 pcs	24 pcs
User Manual	1	1	1

### Preparation of W1 Buffer and Wash Buffer for the first use:


Cat. No:	FABGK000-Maxi (2 preps)	FABGK003 (10 preps)	FABGK003-1 (24 preps)
* Ethanol volume for W1 Buffer	2.5 ml	12 ml	32 ml
**Ethanol volume for Wash Buffer	12 ml	80 ml	160 ml

### Specification:

Principle: spin column (silica membrane)  
 Sample Size: up to 10 ml of fresh/frozen blood  
 up to  $1 \times 10^8$  of cultured cells  
 Column Capacity:  $\leq 500$  µg DNA/column  
 Average DNA yield: 35 µg/ml of whole blood  
 Handling Time: 1 hr  
 Elution Volume: 0.5~1 ml

### Required material to be provided by user

Pipettors and pipet tips  
 Centrifuge: should be capable up to 4,000 xg  
 Thermal incubator  
 Oven (optional)  
 Ethanol (96~100%)  
 Vortex

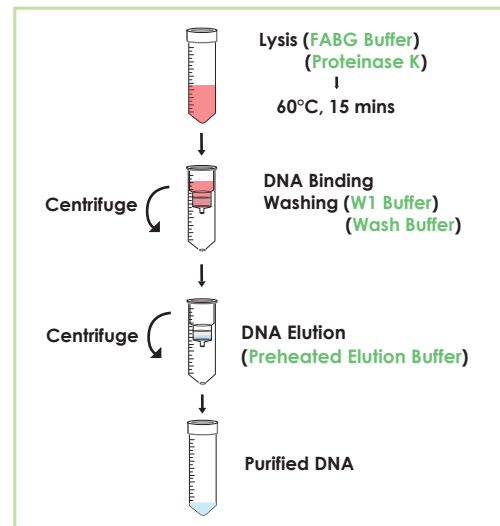

### Important Notes:

1. Buffers provided in this system contain irritants. Wear gloves and lab coat when handling these buffers.
2. Preheat a thermal incubator to 60°C before the operation.
3. Use a centrifuge with a swinging bucket rotor and a force of 4,000~6,000 xg for in all centrifugation steps.
4. Preheat the Elution Buffer or ddH<sub>2</sub>O for step 11 (Elution step).

### Protocol: For Blood DNA Extraction

#### Please Read Important Notes Before Starting the Following Steps.

1. Transfer up to 10 ml sample (whole blood, buffy coat) to a 50 ml centrifuge tube (not provided). If the sample volume is less than 10 ml, add PBS to adjust volume to 10 ml.
2. Add 1000 µl of Proteinase K to the sample and mix well by vortexing. Add 10 ml of FABG Buffer to the sample mixture. **Mix thoroughly by pulse-vortexing.**  
 -Do not add Proteinase K directly to FABG Buffer.
3. Incubate the sample mixture at 60°C for 15 mins to lyse the sample. During incubation, invert the tube every 3~5 mins.
4. **(Optional):** If RNA-free genomic DNA is required, add 80 µl of 100 mg/ml RNase A (not provided) to the sample mixture and incubate at room temperature for 10 mins.
5. Add 10 ml of ethanol (96~100%) to the sample mixture. **Mix thoroughly by vortexing; if precipitate appears, break it by pipetting.**
6. Place a FABG Maxi Column to a 50 ml centrifuge tube (not provided). Transfer 15.5 ml of sample mixture (ethanol added; Including any precipitate) carefully to the FABG Maxi Column. Close the cap and **centrifuge at 4,000~6,000 xg for 3 mins.**

7. Discard the flow-through and transfer the rest sample mixture to the same FABG Maxi Column. Close the cap and **centrifuge at 4,000~6,000 xg for 3 mins** and discard the flow-through.
8. Add 4 ml of W1 Buffer (ethanol added) to the FABG Maxi Column. Close the cap and **centrifuge at 4,000~6,000 xg for 3 mins**. Discard the flow-through and place the FABG Maxi Column back in the 50 ml centrifuge tube.  
-Make sure that ethanol has been added into W1 Buffer at the first open.
9. Add 7 ml of Wash Buffer (ethanol added) to the FABG Maxi Column. Close the cap and **centrifuge at 4,000~6,000 xg for 15 mins**. Discard the flow-through and place the FABG Maxi Column back in the 50 ml centrifuge tube.  
-Make sure that ethanol has been added into Wash Buffer at the first open.  
**-Important Step!** Make sure the residual liquid have been removed completely after centrifugation. It might be necessary to do a further drying by placing the column in a vacuum oven at 70°C for 3 mins.
10. Place the FABG Maxi Column into a new 50 ml centrifuge tube (Elution Tube, provided).
11. Add 0.5~1 ml of preheat Elution Buffer or ddH<sub>2</sub>O (pH 7.5~9.0) to the membrane center of the FABG Maxi Column. **Stand the FABG Maxi Column for 5 mins at room temperature.**  
**-Important Step!** For effective elution, stand the FABG Maxi Column for 5 mins is required to make sure Elution Buffer is absorbed completely by column membrane.
12. **Centrifuge at 4,000 xg for 2 mins** to elute total DNA.  
-Standard volume for elution is 1 ml. If higher DNA yield is required, repeat the DNA Elution step (step 11 and 12) to increase DNA recovery.

### Protocol: For Cultured Cell DNA Extraction

1. Transfer up to  $1 \times 10^8$  of cells to a 50 ml centrifuge tube (not provided). **Centrifuge at 4,000~6,000 xg for 5 mins** to pellet the cells. (If using adherent cells, trypsinize the cells before harvesting.)
2. Resuspend the cells with 10 ml of PBS.
3. Follow the Blood protocol starting from step 2.

### Troubleshooting

Possible reasons	Solutions
<b>Low or no yield of genomic DNA</b>	
<b>Poor cell lysis</b>	
Poor cell lysis because of insufficient Proteinase K activity	Repeat the extraction procedure with a new sample. Make sure the reactive temperature and time is correct.
Poor cell lysis because of insufficient mixing with FABG buffer	Repeat the extraction procedure with a new sample. Mix the sample and FABG Buffer immediately and thoroughly by pulse-vortexing.
Poor cell lysis because of insufficient incubation time	Repeat the extraction procedure with a new sample. Extend the incubation time and make sure that no residual particulates remain.
Ethanol is not added into the lysate before transferring into FABG Maxi Column	Repeat the extraction procedure with a new sample.
<b>Incorrect preparation of Wash Buffer</b>	
Ethanol is not added into W1 and Wash Buffer at the first open	Make sure that the correct volume of ethanol (96~100%) is added into W1 and Wash Buffer at the first open. Repeat the extraction procedure with a new sample.
The volume or the percentage of ethanol is not correct before adding into W1 and Wash Buffer	Make sure that the correct volume of ethanol (96~100%) is added into W1 and Wash Buffer at the first open. Repeat the extraction procedure with a new sample.

Possible reasons	Solutions
<b>Elution of genomic DNA is not efficient</b>	
pH of water (ddH <sub>2</sub> O) for elution is acidic	Make sure the pH of ddH <sub>2</sub> O is between 7.5~9.0.
	Use Elution Buffer (provided) for elution.
Elution Buffer or ddH <sub>2</sub> O is not completely absorbed by column membrane	After Elution Buffer or ddH <sub>2</sub> O is added, stand the FABG Maxi Column for 5 mins before centrifugation.
<b>Column is clogged</b>	
Blood sample contains clots	Repeat the extraction procedure with a new sample. Mix the blood sample well with anti-coagulant to prevent formation of blood clots.
Sample is too viscous	Reduce the sample volume.
<b>Degradation of eluted DNA</b>	
Sample is old	Always use fresh or well-conserved sample for genomic DNA extraction.
Buffer for gel electrophoresis contaminated with DNase	Use fresh running buffer for gel electrophoresis.